

DATASHEET

AmbuMan® Basic is designed to fulfill all requirements for the BLS algorithm training.

And now with the new iQF App, make the training to be given with a high level of precision of adults' CPR performance in cardiac arrest.

CPR becomes the most important factor when a person suffers from cardiac arrest


Brain damage can be caused if no action is taken within four minutes without oxygen under this circumstance.

Therefore, it is important that students are trained in this discipline.

iQF, Training that Works for Life

Helping health instructors to train. The iQF application includes data-driven information and accurate measurements of relevant CPR parameters.

Help your students perform perfect CPR while enjoying games and competing with each other.

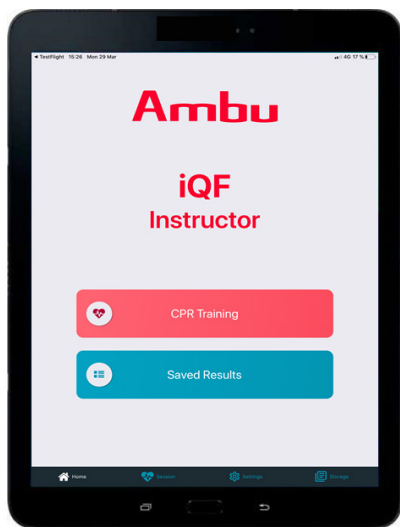

Ambu iQF Instructor

Realistic Ventilation and Correct Compressions

Students will feel how a realistic airway adds resistance, view the rise and fall of the chest when the head is tilted correctly and air is blown into the manikin.

The AmbuMan Basic lets students look, listen and feel the ventilation like in real life.

The unique hygienic system prevents cross contamination.

The correct compression depth and required ventilation volume can be monitored by the feedback indicator on the manikin.

AmbuMan® Basic iQF



KEY BENEFITS

- Robust and reliable cost efficient BLS training manikin with anatomically correct closed Torso
- Visual feedback built into the manikin torso, provides instant feedback of compression amplitude
- Ambu hygienic system for mouth-to-mouth and mouth-to-nose ventilation
- Realistic lifting of the chest during ventilation.
- Easy AED Training

TRAIN IQF, TRAIN FOR LIFE

Visualization in detail of the compression depth and the correct ventilation volume in your fingertips with Instructor and Trainee iQF App. Connected by Bluetooth®.

Ambu

SPECIFICATIONS

Weight Torso	6 kg
Length	70 cm
Width	37 cm
Weight 4 x Torso and Bag	30 kg
Max. Ventilation volume	1200 ml
Max. compression depth	7 cm

MATERIALS

Basic unit	ABS (Acrylonitrile Butadiene Styrene)
Head	Polyvinyl chloride, hard
Chest skin	Polyvinyl chloride, soft
Face piece	Polyvinyl chloride, soft
Headbag disposable	Polyethylene
Trainingsuit	Polyester
Carrying bag	PVC Coated Nylon

ORDERING SPECIFICATIONS

Item no.	
A268008000	AmbuMan Basic iQF incl. 1 AmbuMan Basic Torso with training jacket, carry case, 1 face piece and 25 head bags
A268009000	AmbuMan Basic iQF 4 Pack iQF incl. 4 AmbuMan Basic Torso with training jacket, carry case, 4 face pieces and 100 head bags
234 000 702	Ambu head bags, disposable (set of 100)
234 000 703	Ambu face piece (Set of 5)
234 000 725	Ambu Guedel face piece


* Bluetooth is a registered trademark of Bluetooth SIG, inc.



Compression Indicator



Unique Ambu Hygienic System

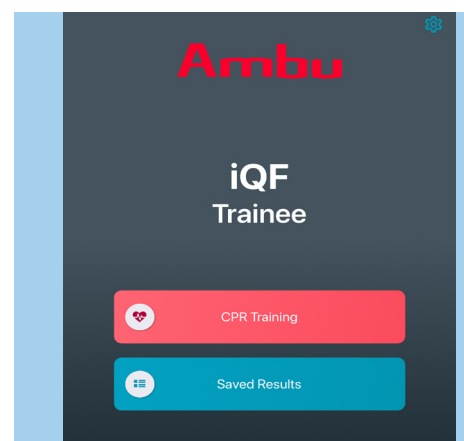

Ambu iQF Trainee

Ambu

■ Ambu A/S
Baltorpbakken 13
2750 Ballerup
Denmark
T +45 72 25 20 00
F +45 72 25 20 50
ambu.com